




BAÑO PARA EXTRACCIÓN DE LÍPIDOS Y PRUEBAS COD


DAIHAN® Extraction Water Bath "WEB", Multi-purpose, 4-6-8 Place, with Certi. & Traceability
With Concentric Ring Cover Sets & Holding Frame, Digital Fuzzy Control, Back Light LCD, up to 100°C, ±0.2°C

CE up to 100°C

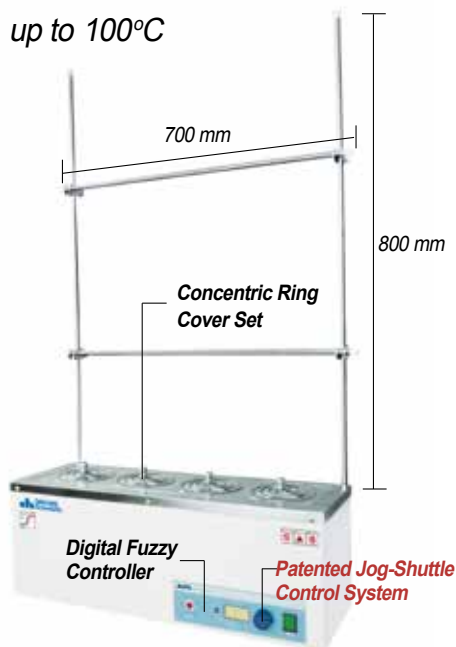

4-Place Extraction Water Bath "WEB-4"

CE up to 100°C

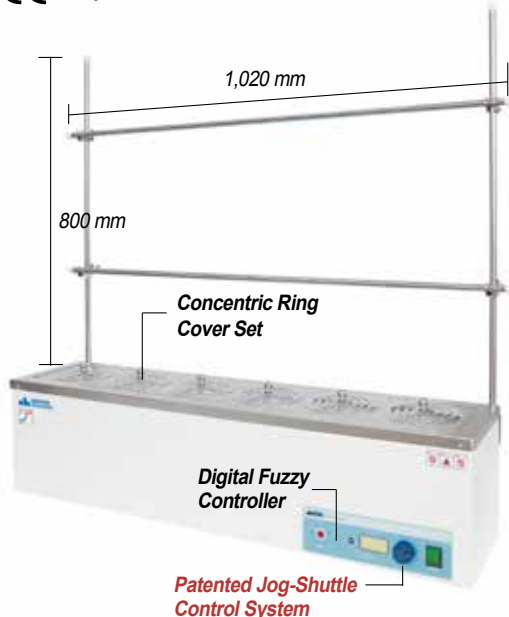

6-Place Extraction Water Bath "WEB-6"



(1) Lid with Concentric Ring Cover Set (Included)

- Made in Korea ■ CE Certified ■ PL (Product Liability) Insurance
 - with Certificate & Traceability : Controlled by Serial Number, Certificate, Delivery-information, and Traceable Data Base System
 - Ideal for Gentle Extraction, Soxhlet Extraction and Determination of COD
 - Digital Fuzzy Control System Implementing Superior Temperature Accuracy
 - Concentric Ring Cover Sets make Effective Steaming Work
 - Patented Jog-Shuttle Control System and Digital LCD with Back-Light Function
- * 4- 6- 8 Extraction Place are available.
 - * Included Concentric Ring Cover Set & Holding Frame
 - * User's Self-compensation Function to Control Difference between Real Temperature and Indicated Temperature (± 10.0°C).
 - * High Quality 304 Stainless-steel Bath and Holding Frame for Superior Durability.
 - * Storage Function for the Set Values of Temperature and Timer.
 - * Alarm Function : Error Status and Timer-end.
 - * Locking Mode Supported for Experimental Safety (Input to Jog-Shuttle can be Disabled).

Specifications & Ordering information

Spec.	Model	"WEB-4" 4 Place	"WEB-6" 6 Place	"WEB-8" 8 Place
Temp. Range & Accuracy		Room Temp. + 5°C to 100°C, ± 0.2°C		
Temp. Uniformity & Sensor		± 1.0°C, PT100		
No. of Extraction Holes		4 Place	6 Place	8(4x 2) Place
Heating Power		1.4 kW	1.8 kW	1.4 kWx 2
Timer & Alarm		99hr 59min (Delay & Continuous run), Error status & Timer-end		
Display		Digital LCD with Back-Light Function		
Controller		Digital Fuzzy Controller with Jog-Shuttle Switch (Turn+Push)		
Safety Circuit		Over-Temp./-Current protection, Sensor Error Detection, Leakage Breaker		
Material		Stainless-steel (#304)-Internal, Powder Coated steel-External, Stainless-steel Holding Frame		
Flask Capacity		250 ~ 500ml Flask		
Lid		Concentric Ring Cover Set Included		
Top Opening Dia.		Φ 36 Φ 65 Φ 96 Φ 110mm		
Bath Depth		122mm		
Drain Valve		Mounted		
Dimension (w×d×h)	Internal	610x 160x h160 mm	910x 140x h160 mm	610x 290x h160 mm
	External	680x 245x h275 mm	975x 245x h275 mm	720x 375x h275 mm
Packing Size (mm) and Gross Weight (kg)		1,100x 350x h400, 18 kg	1,250x 350x h400, 25 kg	900x 450x h400, 34 kg
Electrical Reqtts.		1 Phase, AC 120V, 60Hz or 230V, 50 / 60Hz		
Cat. No.	230V	DH.WEB01004	DH.WEB01006	DH.WEB01008
	120V	DH.WEB02004	DH.WEB02006	DH.WEB02008

* Other Specifications are available upon Customer's Request.

Spare Lid with Concentric Ring Cover Set and Holding Frame Set

- DH.WEB11000 (1) Lid with Concentric Ring Cover Set, "EBC100", for all "WEB" models
- DH.WEB21004 Holding Frame Set, "EBF100", for "WEB-4", 2eax 800mm Pole & 2eax 700mm Pole
- DH.WEB21006 Holding Frame Set, "EBF200", for "WEB-6", 2eax 800mm Pole & 2eax 1,020mm Pole
- DH.WEB21008 Holding Frame Set, "EBF300", for "WEB-8", 4eax 800mm Pole & 4eax 740mm Pole

