



## STOMACHER 400ml

**BagMixer®** 400 P is an easy-to-use and powerful lab blender. Adapted to all kinds of applications and with a guarantee of optimal bacterial extraction, it is a great tool for lab analyses.



### TECHNICAL SPECIFICATIONS

- All stainless steel 304L
- Plain door
- Useful volume: 50 - 400 mL
- Fixed speed (8 strokes/s)
- Variable blending time (30 - 210 s or ∞)
- Q-Tight® closing: tight closing for no more leaks from the top of the bag
- Fully opening door (270°)
- Also available in CoolMix® version
- Voltage and frequency: 110-127V~ / 220-240V~ 50-60Hz
- Fuses: T 3.15A L 240V / T 1.25A L 5x20 mm
- Max power: 135 W
- Dimensions (w x d x h): 27 x 38.5 x 28 cm, weight: 16.6 kg
- Box (w x d x h): 49.5 x 35 x 36.5 cm, weight: 18.6 kg
- In compliance with: ISO 7218, ISO 6887-1 and FDA BAM (Bacteriological Analytical Manual)
- Manufactured under CE, RoHS, WEEE - Products made for INTERSCIENCE by Interlab, an ISO 9001 certified company
- Designed and made in France

### Best-seller

- Sold in more than 90 countries
- 35 years experience

### High-quality

- Optimal analysis with the sterile **BagFilter®** and **BagPage®**
- In compliance with ISO 7218 standard
- Made in France by interscience

Standards:



Manufactured under:



Product made for INTERSCIENCE by Interlab, an ISO 9001 certified company

### Brushless motor

- No wearing parts
- - 50% power consumption
- Extra silent

### Safety

- The handle closes the blender and starts the blending
- **Q-Tight®**: no leaks from the top of the bag
- Security drip tray

PATENTED

### Click & Clean® removable paddles

Remove the paddles and easily clean the blending chamber



### Easy cleaning

- All stainless steel
- Door opens at 270°
- Removable **Click & Clean®** paddles

### Guarantee

- Shock afterlife time and window door guarantee
- 3 year warranty (after registration)

PATENTED


### Adjustable blending power

- Up to 28 kg of pressure: the hardest samples don't resist! e.g.: nuts, pills, seeds, dog biscuits etc.



### Optimal extraction

- From as little as 30 seconds
- Easy control with the digital screen (speed, time, power)

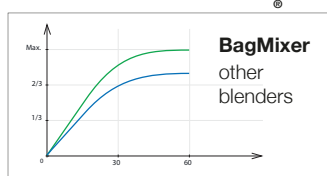

### OPTIONS / ACCESSORIES

- Ref 026 000 Drip-tray: Stainless steel security drip tray for BagMixer® 400

### DELIVERED WITH

- 1 power cable
- 1 user's manual
- 1 quick user guide
- 1 free pack of blender bags





The **BagMixer**<sup>®</sup> lab blender guarantees optimal extraction by the action of its two paddles. Adapted to any sample (food, agricultural, chemical, medical, environmental) to be analyzed between 50 mL and 400 mL, it guarantees homogenization in 30 to 60 seconds.



**BagMixer**<sup>®</sup> 400 P  
> The essential



**BagMixer**<sup>®</sup> 400 W  
> Window door



**BagMixer**<sup>®</sup> 400 CC  
> Click & Clean<sup>®</sup>



**BagMixer**<sup>®</sup> 400 S  
> Silent blending



**BagMixer**<sup>®</sup> 400 SW  
> Window door & LiquidSensor

Also available: packs including bags and accessories; 100 mL and 3500 mL lab blenders. For more info, please contact us.

For more info on **BagMixer**<sup>®</sup> S & SW models, please refer to **BagMixer**<sup>®</sup> S & SW brochure.

## BagMixer<sup>®</sup> packs:

### Starter Pack



- 1 **BagMixer**<sup>®</sup>
- + **Sterile bags:** 500 BagFilter<sup>®</sup>  
+ 500 BagLight<sup>®</sup>
- + **Accessories:** 50 BagClip<sup>®</sup>  
+ 1 BagOpen<sup>®</sup> + 1 BagRack<sup>®</sup>

### Classic Pack




- 1 **BagMixer**<sup>®</sup>
- + **Sterile bags:** 500 BagFilter<sup>®</sup>  
+ 500 BagLight<sup>®</sup> + 500 BagPage<sup>®</sup>
- + **Accessories:** 50 BagClip<sup>®</sup>  
+ 1 BagOpen<sup>®</sup> + 1 BagRack<sup>®</sup>  
+ 1 PetriPile<sup>®</sup> + 1 BagPipet<sup>®</sup>  
+ 1 000 BagTips<sup>®</sup> Regular  
+ 6 Clinet<sup>®</sup> 101

### Premium Pack




- 1 **BagMixer**<sup>®</sup> & 1 DiluFlow<sup>®</sup>
- + **Sterile bags:** 1 000 BagFilter<sup>®</sup>  
+ 2 000 BagLight<sup>®</sup> + 1 000 BagPage<sup>®</sup>
- + **Accessories:** 50 BagClip<sup>®</sup>  
+ 1 BagOpen<sup>®</sup> + 2 BagRack<sup>®</sup>  
+ 2 PetriPile<sup>®</sup> + 1 BagPipet<sup>®</sup>  
+ 1 000 BagTips<sup>®</sup> Regular  
+ 6 Clinet<sup>®</sup> 101

## BagSystem<sup>®</sup> complete range for sample preparation before analysis

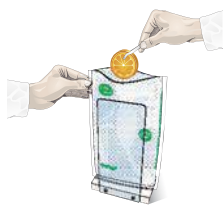

**BagFilter**<sup>®</sup>  
Collect

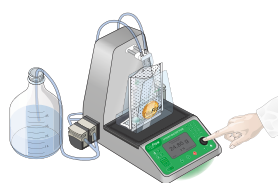

**DiluFlow**<sup>®</sup>  
Dilute



**BagMixer**<sup>®</sup>  
Blend



**BagTools**<sup>®</sup>  
Pipet



# Technical specifications

		BagMixer® 400 P	BagMixer® 400 W	BagMixer® 400 CC		
<b>Key features</b>	Removable and autoclavable paddles, <b>Click &amp; Clean®</b> system	-	-			
	Blending power indicator	-	-			
	Adjustable blending power (up to 28 kg of pressure)	-	-			
	Digital screen	-	-			
	Leak free <b>Q-tight®</b> system	-	-			
	Side-by-side paddle stop	-	-			
	Silent brushless motor	-	-			
	Window door	-	-			
	Security drip tray	in option				
	270° opening door					
	All stainless steel					
	Compact and ergonomic					
<b>Technical specifications</b>	Door weight	1050 g	950 g			
	Blending time	30, 60, 90, 120, 150, 180, 210 sec., ∞		from 1 sec. to 59 min, ∞		
	Blending speed	fixed speed: 8 strokes/sec.		variable speed: 4, 6, 8, 10 strokes/sec.		
	Dimensions (w x d x h)	27 x 38.5 x 28 cm		27 x 41.5 x 26.5 cm		
	Weight	16.6 kg		16 kg		
	Power	110-127/220-240 V~ 50-60Hz		100-240 V~ 50-60Hz		
	Capacity	50-400 mL				
<b>Warranties</b>	Made in France					
	Life-time warranty	Unbreakable shock absorbers and window door				
	3 year warranty (after registration)					
<b>Manufactured in compliance with</b>						

## Discover our complete range for microbiology



**BagFilter®**  
Sterile bags



**DiluFlow®**  
Gravimetric dilutors



**FlexiPump®**  
Peristaltic pumps



**easySpiral Dilute®**  
Automatic diluter et platers



**Scan®**  
Automatic colony counters  
& inhibition zone readers