

BACTI-CINERATOR STERILIZER MAX. 825°C

REF: HY800




Description

HY-800 (small-diameter) Bacti-cinerator Sterilizer uses infrared heat sterilization, because it is easy to use, simple to operate, no open flame, no fear of wind, substitute alcohol lamp, safe to use, can be widely used in the environment of biosafety cabinet, purification workbenches, suction fans, and mobile vehicles, for microbiological experiments.

Features

- Based on the principle of electronic infrared, the products are widely used in the sterilization of inoculation rings or inoculation needles, completely replacing alcohol lamps.
- The temperature inside the heating hole can reach 800°C or above, sterilization takes only 5-7 seconds and sterilization is through. The product can be used in an anaerobic chamber.
- Small volume, light weight, appearance and easy to clean. Heating body can adjust angle, easy operation.
- Organic substances in the ceramic funnel pipe depths ashing prevent infectious spatter and cross contamination.
- Precise temperature control technology, avoid the heating components work overload, long life.

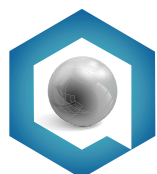

BACTI-CINERATOR STERILIZER MAX. 825°C

REF: HY800

Technical specifications

Parameters	Description
Center Max. Temp	825°C ± 50°C
Low level Temp	480°C
Heating Speed	20 min (to max. Temp.)
Max. Sterilize diameter	Ø14mm
Heating element lenght	140mm
Adjustable angle	75° - 120°
Voltage	AC110V / 60HZ
Poder	230W
Fuse	250V, 2A/3A, Ø5X20
Dimension (WxDxH)	87x151x160mm
Net Weight	1kg