


ULTRASONIC SPARK PLUG CLEANER 6.8L

REF: UC-9180L



Description

- One-key degassing can quickly remove air in water, and the ultrasonic cleaning effect is better
- Independently developed MCU-sweep ultrasonic generator drive, uniform and strong ultrasonic effect
- Independently developed industrial high Q value transducer with high ultrasonic conversion efficiency and long service life
- Power adjustment function, can meet different cleaning requirements
- Imported high strength glue, seedless nail bonding process, higher ultrasonic conversion efficiency
- Digital timing, cleaning time 10s-100min and digital heating control, temperature 20-80°C arbitrary setting
- The inner groove adopts imported SUS304, 1.0mm stainless steel stamping, ultrasonic exchange effect is better, longer service life
- The washing basket is made of high quality 304 stainless steel mesh welding, surface electrolytic polishing treatment
- The shell is made of high quality stainless steel, with better anticorrosion effect

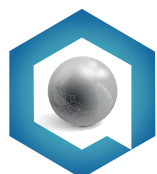

QUIMICOMPANY

EQUIPOS E INSUMOS PARA LABORATORIO

denmark
SCIENTIFIC

ULTRASONIC SPARK PLUG CLEANER 6.8L

REF: UC-9180L

LCD display and control, simple operation



Safety switch, power interface



Portable drain valve



Foot pad, heat dissipation hole design



Stainless steel one-piece inner groove



Plastic handle

☎ (601) 300 3161 + (57) 316 532 6905

🌐 www.quimicompany.com.co

✉ info@quimicompany.com.co

ULTRASONIC SPARK PLUG CLEANER 6.8L

REF: UC-9180L

Technical specifications

Parameters	Description
The tank size(L*W*H) mm	300*150*150
Dimensions(L*W*H) mm	367*176*292
Capacity(L)	6.8
Gross weight (KG)	6.2
Ultrasonic power (W)	180
Heating power(W)	200
Drain valve	Yes
Frequency(KHz)	39±1.5KHz
Time (Min)	10Sec~100Min
The heating temperature (°C)	20~80
Whashing basket	Included 1 Unit.